

CHAIR

OP-S8

EXCELLENCE IN ENT CHAIRS

- Resistant
- Comfortable
- Totally electrical





producing quality

[OP-S8]

FOLDING ARMRESTS

They are composed of a stainless steel axle padded with washable polyurethane. They fold down completely to provide more comfort for patients with reduced mobility and for lateral observations and practices.



STRONG HEADREST

Perfectly designed to keep the head of the patient in the required position during the examination.


Adjustable in height and inclination, it incorporates a quick removal system which facilitates the execution of positional tests and routine practices.



COMFORTABLE FOR DOCTOR AND PATIENT

AN ERGONOMIC PRODUCT

The OPTOMIC chairs design is based on the search for the comfort of the patient and ergonomics during exploration. This is the reason the geometry of its structure adapts to achieve the optimum distance between patient and doctor in every step of examination and treatment.

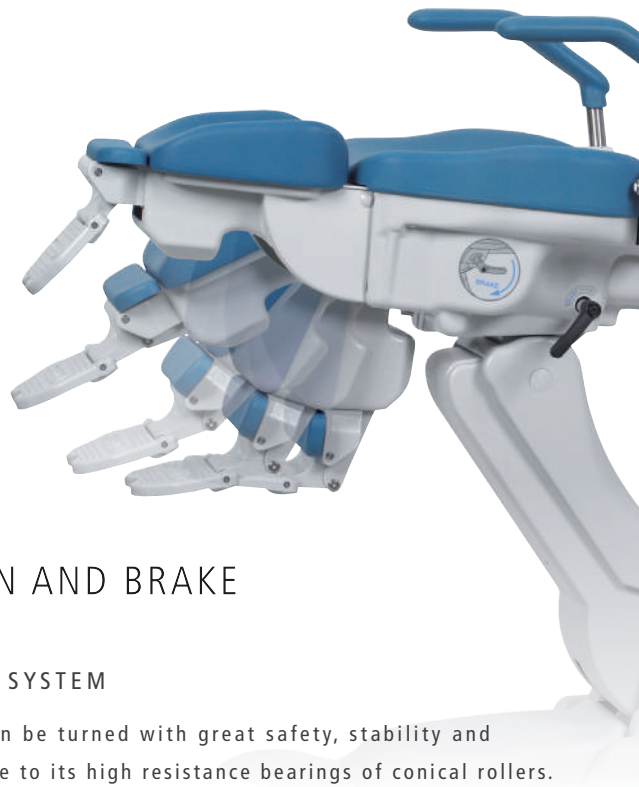

ROTATION AND BRAKE

ACCESSIBLE SYSTEM


The patient can be turned with great safety, stability and smoothness due to its high resistance bearings of conical rollers.

It has a double lever braking system, placed on either side of the seat, to avoid the doctor wasting time moving around the patient.

The 300° turn facilitates the examination of both ears.



CHAIR OP-S8

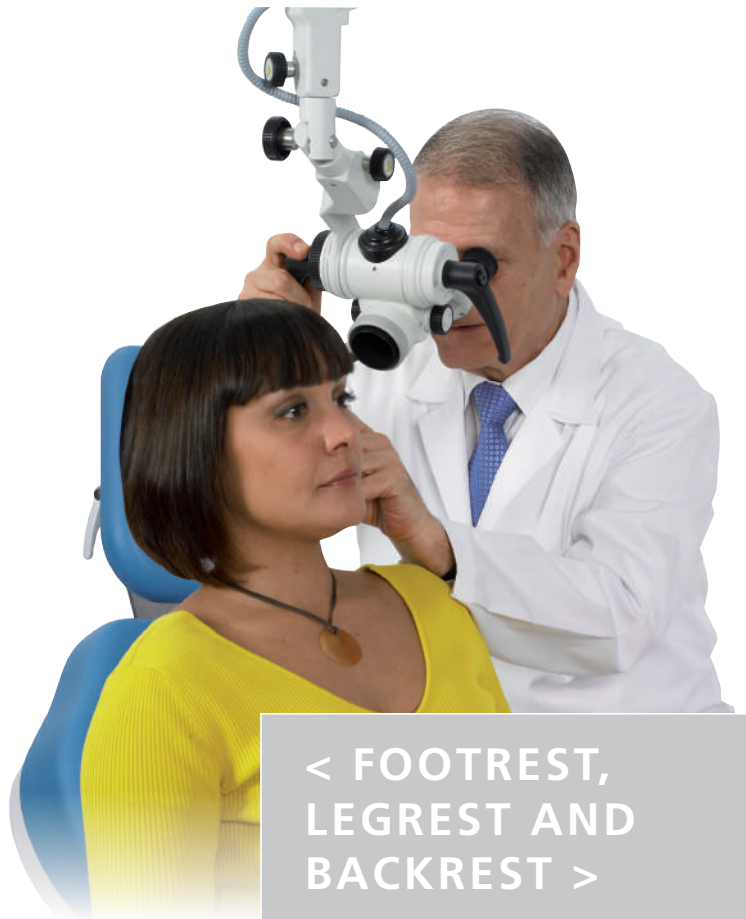

< FOOTREST,
LEGREST AND
BACKREST >

MOTORIZED ELEVATION SYSTEM

A PRECISE SYSTEM

A pedal facilitates the perfect control of the chair. It activates the up /down and inclination movements. It has two memories for user personalization.



HIGH RESISTANCE AND SMALL SIZE

The OP-S8 chair has a completely electrical system for the movement of the footrest and legrest, which is synchronized.

The legrest is concave to hold the patients calves for greater comfort and relaxation.

The range of the backrest movement goes from 80° to the Trendelenburg position.



producing quality

The OPTOMIC OP-S8 chair is designed and manufactured using the best materials and most modern technologies. OPTOMIC's philosophy is based on the idea that both patient and doctor must be able to find a comfortable and safe position, to avoid unnecessary fatigue during an examination or surgical treatment.

CHAIR OP-S8



OP-S8



UP / DOWN MOVEMENT
 BACKREST AND LEGREST MOVEMENT
 PROGRAMMABLE MEMORIES
 MOVEMENT CONTROL
 UPHOLSTERY
 FIRE PROOF CLASS
 FOOTREST
 ARMRESTS
 HEADREST

CHAIR ROTATION
 TRENDENBURG POSITION
 BACKREST ANGLE
 SEATING HEIGHT
 MAINS CONNECTION
 DIMENSIONS (BETWEEN BASE TRIMS)
 MAXIMUM LOAD

Electrical
Electrical
7 memories (6 programmable + Trendelenburg)
2 keypads on the backrest + Foot-operated joystick
PUR melting, no seams, latex free, washable
Fire proof N°UL94HB
Metallic with automatic fold-down
Folding armrests
Adjustable headrest, variable height and inclination. Removable.
300° with brake on either side
-10° inclination
From 80° to 190°
Max. 900 mm - Min. 520 mm
110/240v - 50-60Hz
600mm (w) x 730mm (h)
200Kg

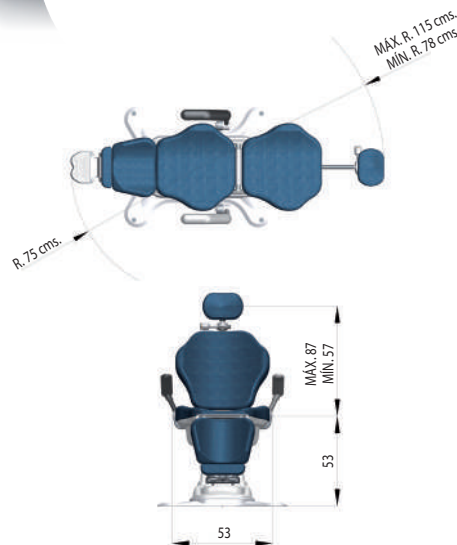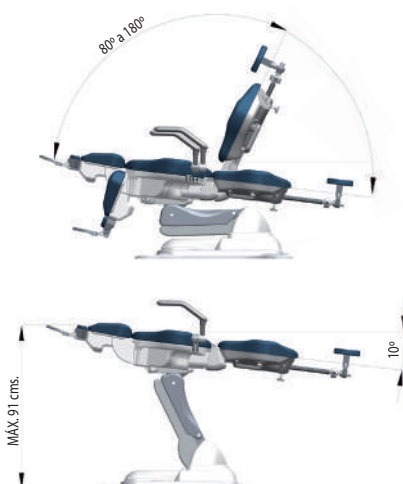

OPTOMIC ESPAÑA S.A.

Madroño 4. E-28770, Colmenar Viejo. Madrid - Spain
 Tel: +34 91 804 44 44
 export@optomic.com
 www.optomic.com



Specifications are subject to change without notice

© Optomic - Ref. C082ENG22

